

**Conditions:**

Injection: Split 250°C, 0.5µL, 1:50 Split Ratio.

Sample dilution: 1mg/1mL each in n-Hexane.


**Contact us** for additional information.

**Menthols**

Chiral Separation on MEGA-DEX Columns

Column: **MEGA-DEX DMT Beta - 0.25mm, 0.25µm, 25m**

Catalog Code: DEX-DMT-B-025-025-25

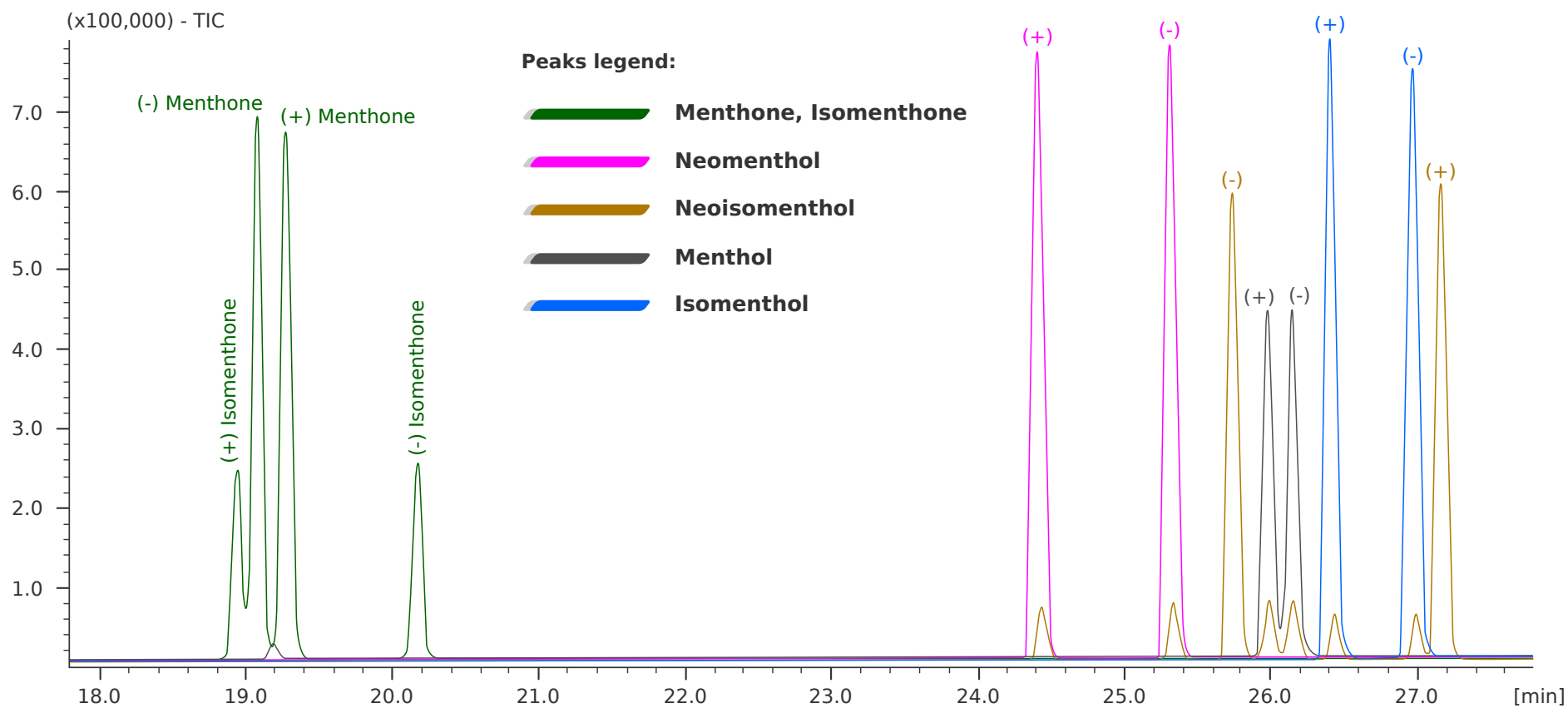

Acknowledgement: Prof. C.Bicchi, C.Cordero et al. - Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via Pietro Giuria, 9 - Torino - Italy.

