

Pesticides

Column: **MEGA-1701 FAST - 0.10mm, 0.10µm, 5m**

Catalog Code: F-1701-010-010-5

Conditions:

Injection: Split 230°C, 0.5µL, 1:200 Split Ratio.

Sample Dilution: 1% in Cyclohexane.


Detector: FID 250°C.

Oven Temp.: 50°C (0.1min), 40°C/min, 250°C.

Carrier Gas: Hydrogen, 0.5 mL/min.

Peaks:

1. Alpha-HCH
2. Gamma-HCH
3. Heptachlor
4. Chlorothalonil
5. /
6. Methyl Parathion
7. Malathion
8. Fenotrothion
9. Ethyl Parathion
10. /
11. Fenitrothion
12. Chlordane-Cis + Trans
13. Dieldrin
14. o,p'-DDT
15. Beta-Endosulfan
16. p,p'-DDT
17. /
18. Tetradifon

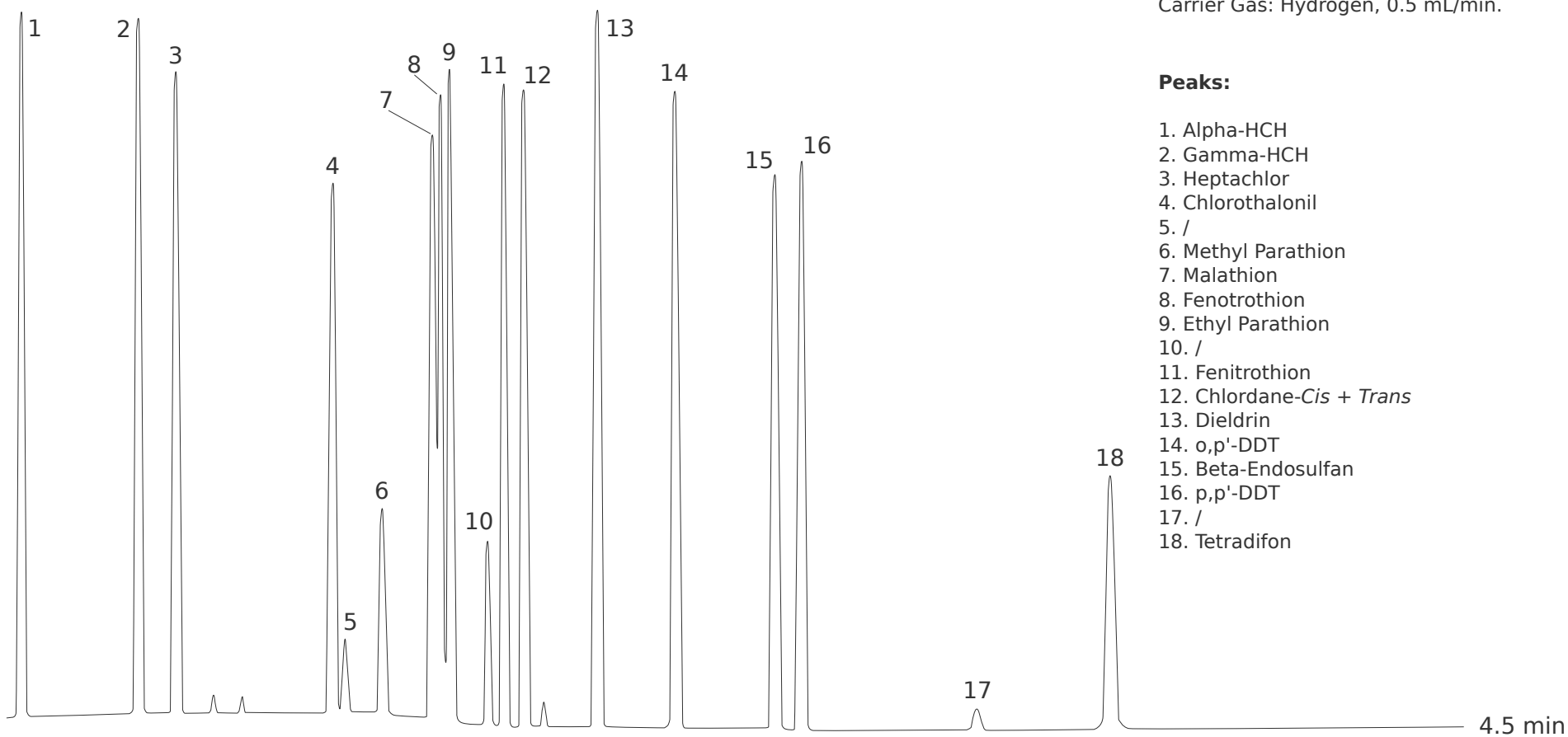

Acknowledgement: Prof. C.Bicchi, C.Brunelli et al. - Università di Torino, Dipartimento di Scienza e Tecnologia del Farmaco, Via Pietro Giuria, 9 - Torino - Italy.

